

ACI 03 ACIDS

Analysis of 15 Organic Acids on BP20

KEYWORDS:

Organic Acids


Fatty Acids

ORGANIC ACIDS

Part No.: 054433
 Phase: BP20, 0.25µm
 Column: 30m x 0.32mm ID
 Initial Temp.: 70°C
 Rate: 10°C/min
 Final Temp.: 260°C, 5 min
 Detector: FID
 Injection Mode: Split
 Carrier Gas: H₂, 6psi

Components

1. Valeric Acid
2. Caproic Acid
3. Heptanoic Acid
4. Octanoic Acid
5. Nonanoic Acid
6. Decanoic Acid
7. Undecanoic Acid
8. Dodecanoic Acid
9. Tridecanoic Acid
10. Tetradecanoic Acid
11. Pentadecanoic Acid
12. Hexadecanoic Acid
13. Heptadecanoic Acid
14. Octadecanoic Acid
15. Arachidic Acid

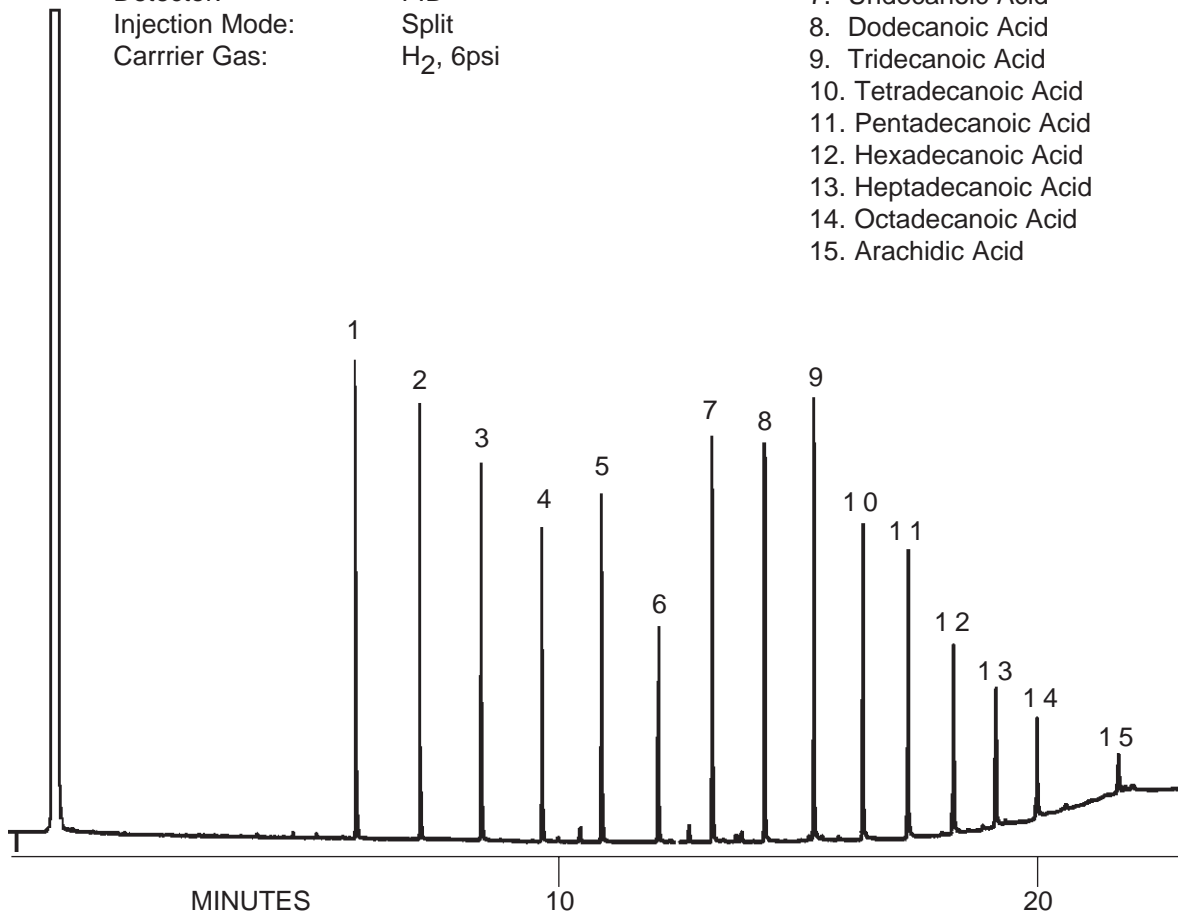

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email: support@sge.com

SGE Incorporated (USA)
email: usa@sge.com

SGE (France) Sarl
email: france@sge.com

SGE Deutschland GmbH
email: germany@sge.com

SGE Japan Inc.
email: sgejpn@itjit.ne.jp

SGE Australia Pty. Ltd.
email: sales@sge.com

SGE Europe Ltd. (UK)
email: uk@sge.com

SGE (Italia) Srl.
email: sge.italia@tin.it

SGE China Service Centre
email: sge@bj.col.com.cn

SGE India
email: sgeindia@bom4.vsnl.net.in