

## Broad Range Pesticide Mixture on BPX5

**KEYWORDS:**


Organochlorine pesticides  
OCP  
GC/MS

**Components**

1. Methyldymron decomp. I
2. Methyldymron decomp. II
3. Trichlorfon (DEP)
4. Etridiazole
5. Chloroneb
6. Benfluralin
7. Pencycuron
8. Simazine (CAT)
9. Diazinon
10. Propyzamide
11. Chlorothalonil (TPN)
12. Terbutcarb (MBPMC)
13. Tolclofos - methyl
14. Fenitrothion (MEP)
15. Chlorpyrifos
16. Pendimethalin
17. Isophenphos
18. Methyldymron
19. Captan
20. Butamifos
21. Napropamide
22. Flutolanil
23. Isoprothiolane
24. Isoxathion
25. Bensulide
26. Mepronil
27. Pyridaphenthion
28. Iprodione

**Broad Range Pesticides**

Part No.: 054142  
Phase: BPX5, 0.25µm  
Column: 30m x 0.22mm  
Initial Temp.: 50°C, 1min  
Rate 1: 30°C/min  
Temp 2: 190°C  
Rate 2: 4°C  
Final Temp.: 280°C, 5min  
Detector: 5973 MSD  
Injector Mode: Splitless, 1.0min  
Carrier Gas: Helium, 35cm/sec  
Injection Volume: 1µL  
Interface/Source/Quad: 320°C, 250°C, 150°C

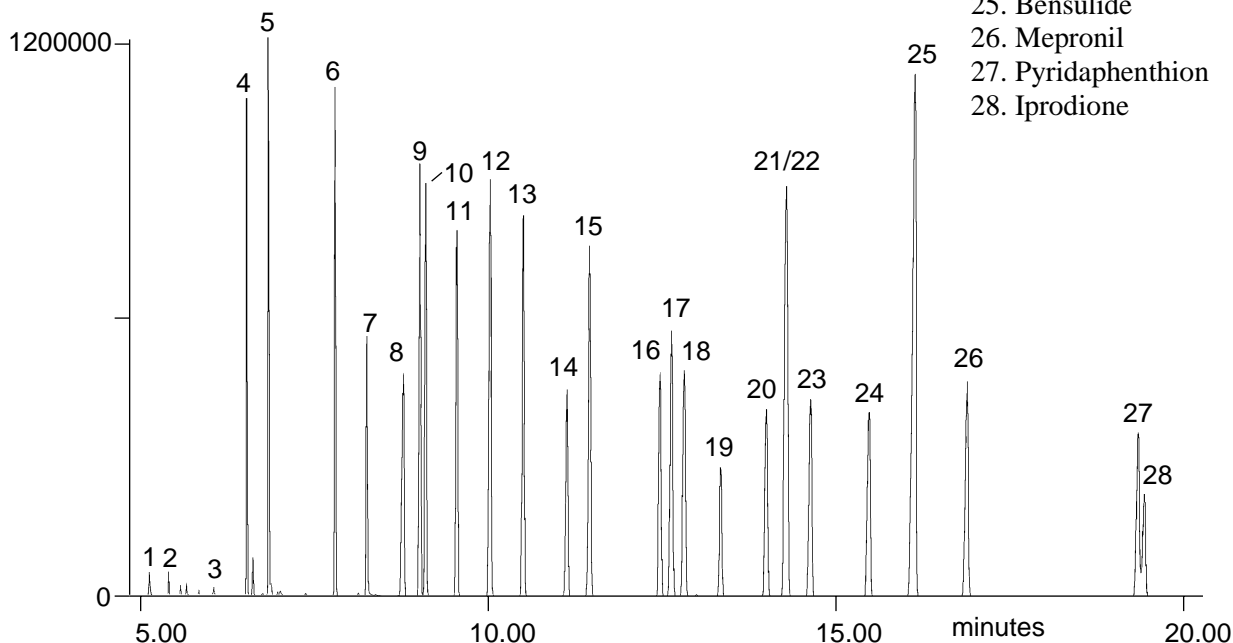

visit [www.sge.com](http://www.sge.com) or contact our technical Customer Support Team on: [techsupport@sge.com](mailto:techsupport@sge.com)