


Some commonly used Herbicides/Pesticides on BPX50

Part No.: 054751
Phase: BPX50, 0.25µm
Column: 30m x 0.25mm ID
Initial Temp.: 50°C, 1min
Rate 1: 70°C/min
Temp 2: 190°C
Rate2: 4°C/min
Final Temp.: 280°C, 5min
Detector: HP5973 MSD
Injector Mode: Splitless, 1.0min
Carrier Gas: Helium, 20psi
Injection Volume: 1µl

Components

1. Methylmymron decomp. I
2. Methylmymron decomp. II
3. Trichlorofon (DEP)
4. Etridiazole
5. Chloroneb
6. Benfluralin
7. Pencycuron
8. Simazine (CAT)
9. Diazinon
10. Propyzamide
11. Chlorothalonil (TPN)
12. Terbucarb (MBPMC)
13. Tolclofos
14. Fenitrothion (MEP)
15. Chlorpyrifos
16. Pendimethalin
17. Isophenphos
18. Methylmymrom
19. Captan
20. Butamifos
21. Napropamide
22. Flotolanil
23. Isoprothiolane
24. Isoxathion
25. Mepronil
26. Pyridaphenthion
27. Iprodione

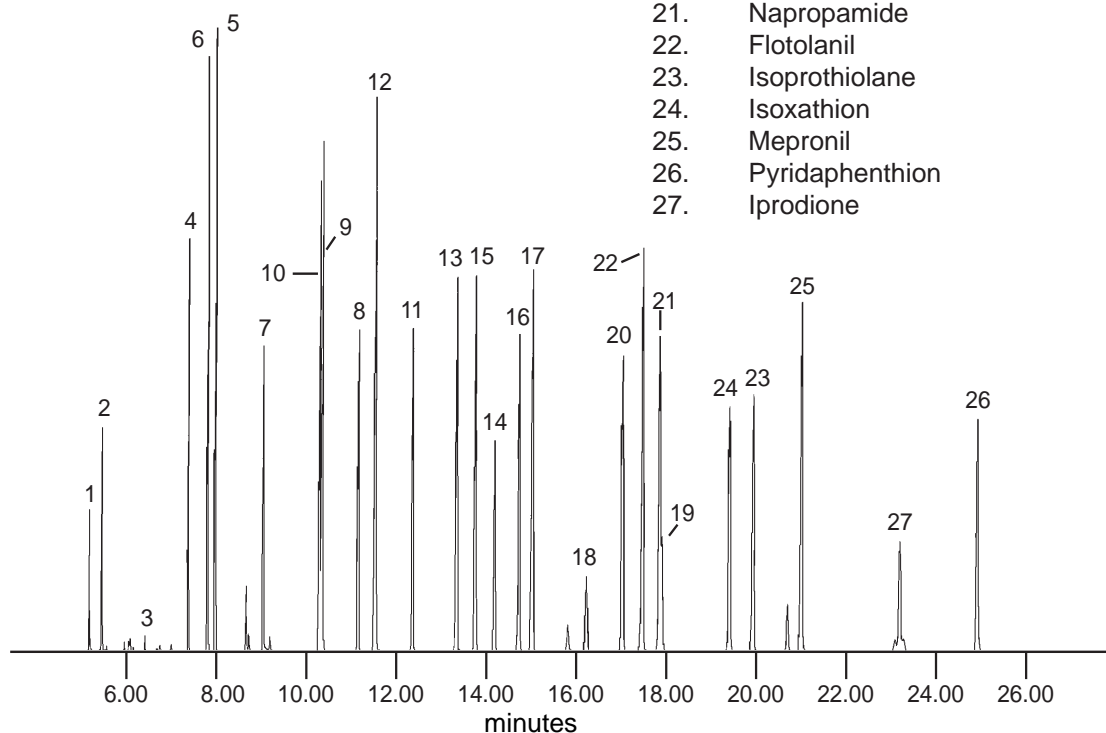

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email: support@sge.com

SGE Incorporated (USA)
email: usa@sge.com

SGE (France) Sarl
email: france@sge.com

SGE Deutschland GmbH
email: germany@sge.com

SGE Japan Inc.
email: sgejpn@itjit.ne.jp

SGE Australia Pty. Ltd.
email: sales@sge.com

SGE Europe Ltd. (UK)
email: uk@sge.com

SGE (Italia) Srl.
email: sge.italia@tin.it

SGE China Service Centre
email: sge@bj.col.com.cn

SGE India
email: sgeindia@bom4.vsnl.net.in