

FOO 02
FOOD


Analysis of Food Additives -Antimicrobials on BP5

FOOD ADDITIVES - ANTIMICROBIALS

Part No.: 054186
Phase: BP5, 0.5µm film
Column: 25m x 0.32mm ID
Initial Temp: 160°C, 0 min
Rate: 15°C/min
Final Temp: 280°C, 0min
Detector: FID
Sensitivity: 256 x 10⁻¹² AFS
Injection Mode: Split

Components

1. Methyl Paraben
2. Ethyl Paraben
3. Propyl Paraben
4. Butyl Paraben
5. Heptyl Paraben

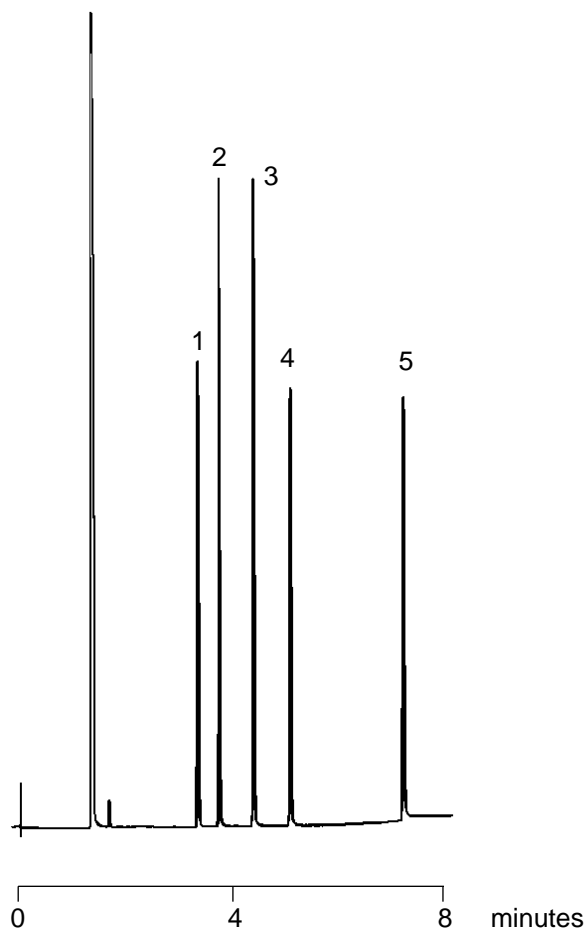

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email: support@sge.com

SGE Incorporated (USA)
email: usa@sge.com

SGE (France) Sarl
email: france@sge.com

SGE Deutschland GmbH
email: germany@sge.com

SGE Japan Inc.
email: sgejpn@itjit.ne.jp

SGE Australia Pty. Ltd.
email: sales@sge.com

SGE Europe Ltd. (UK)
email: uk@sge.com

SGE (Italia) Srl.
email: sge.italia@tin.it

SGE China Service Centre
email: sge@bj.col.com.cn

SGE India
email: sgeindia@bom4.vsnl.net.in