

FOO 10

FOOD

Standard mixture of Triglycerides on BPX5

KEYWORDS:

Triglycerides

Notes: This analysis requires on-column injection to provide quantitative results


Standard Triglycerides

STANDARD MIXTURE OF TRIGLYCERIDES

Part No.: 054133
 Phase: BPX5, 0.25µm
 Column: 12m x 0.53mm I.D.
 Initial Temp.: 100°C, 0.5 min.
 Rate 1: 50°C/min.
 Temp.: 280°C
 Rate 2: 10°C/min.
 Final Temp.: 360°C, 5 min.
 Carrier Gas: Helium, 6psi
 Injection Mode: On-column (OCI-5)
 Injection Volume: 1µl

Components

| PEAK No. | TRIGLYCERIDE CARBON No. |
|----------|-------------------------|
| 1 | T36 |
| 2 | T38 |
| 3 | T40 |
| 4 | T42 |
| 5 | T44 |
| 6 | T46 |
| 7 | T48 |
| 8 | T50 |
| 9 | T52 |
| 10 | T54 |

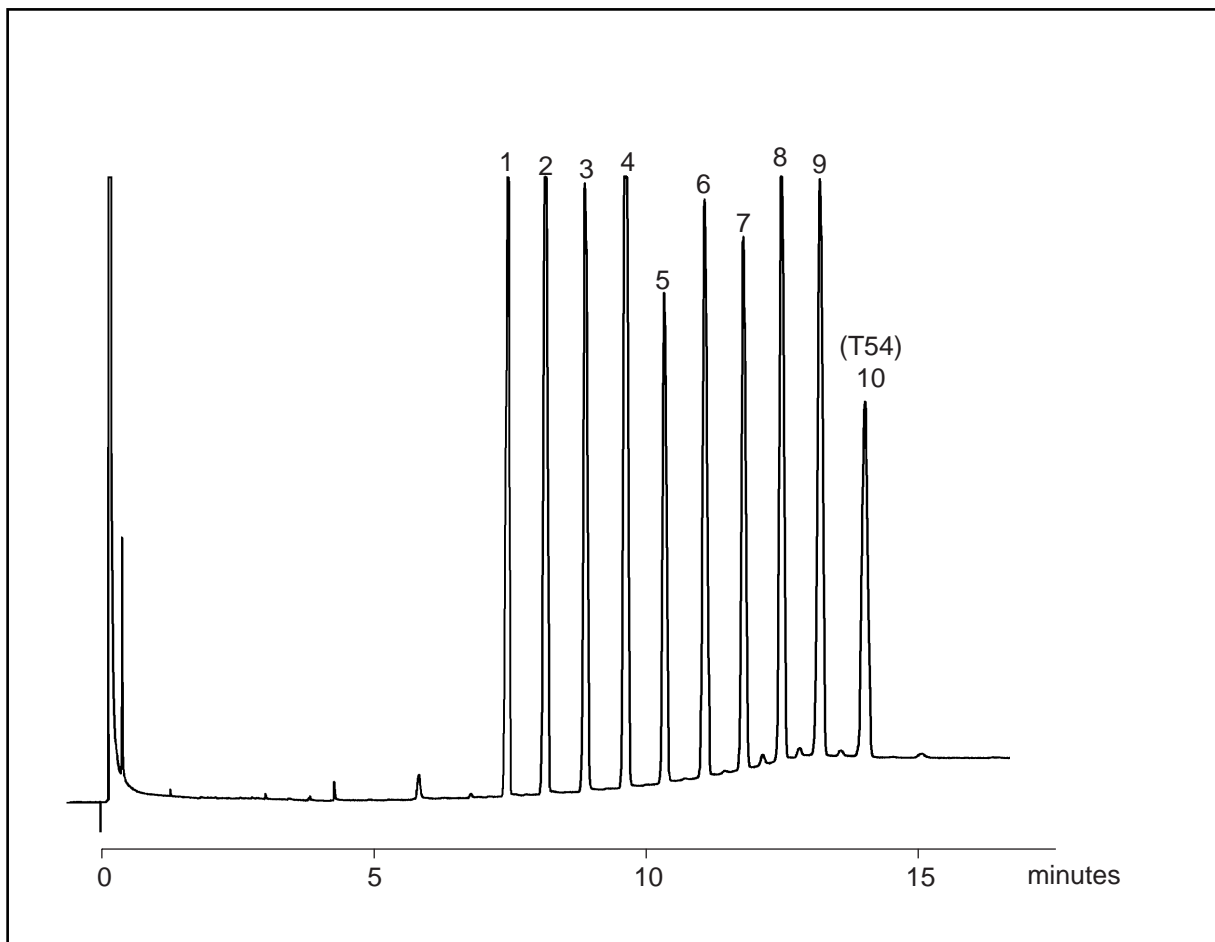

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email: support@sge.com

SGE Incorporated (USA)
email: usa@sge.com

SGE (France) Sarl
email: france@sge.com

SGE Deutschland GmbH
email: germany@sge.com

SGE Japan Inc.
email: sgejpn@tjit.ne.jp

SGE Australia Pty. Ltd.
email: sales@sge.com

SGE Europe Ltd. (UK)
email: uk@sge.com

SGE (Italia) Srl.
email: sge.italia@tin.it

SGE China Service Centre
email: sge@bj.col.com.cn

SGE India
email: sgeindia@bom4.vsnl.net.in