

FOO 22
FOOD

Cod Oil Fatty Acids (Underivatised) on BP21


KEYWORDS:

Cod Oil

Free Fatty Acids

Part No.: 054468
Phase: BP21, 0.25µm
Column: 25m x 0.32mm ID
Initial Temp.: 70°C, 0.5min
Rate 1: 6°C/min
Final Temp.: 260°C, 10min
Detector: FID, 280°C
Injector Mode: Split, 50:1
Carrier Gas: Helium, 15psi
Injection Volume: 1µL
Injection Temp: 280°C

- Components**
- 1. C14:0
 - 2. C16:0
 - 3. C16:1
 - 4. C18:0
 - 5. C18:1
 - 6. C18:2
 - 7. C18:3
 - 8. C20:0
 - 9. C22:0
 - 10. C24:0
 - Int. 1. 3-Methylbutanoic Acid
 - Int. 2. Octanoic Acid

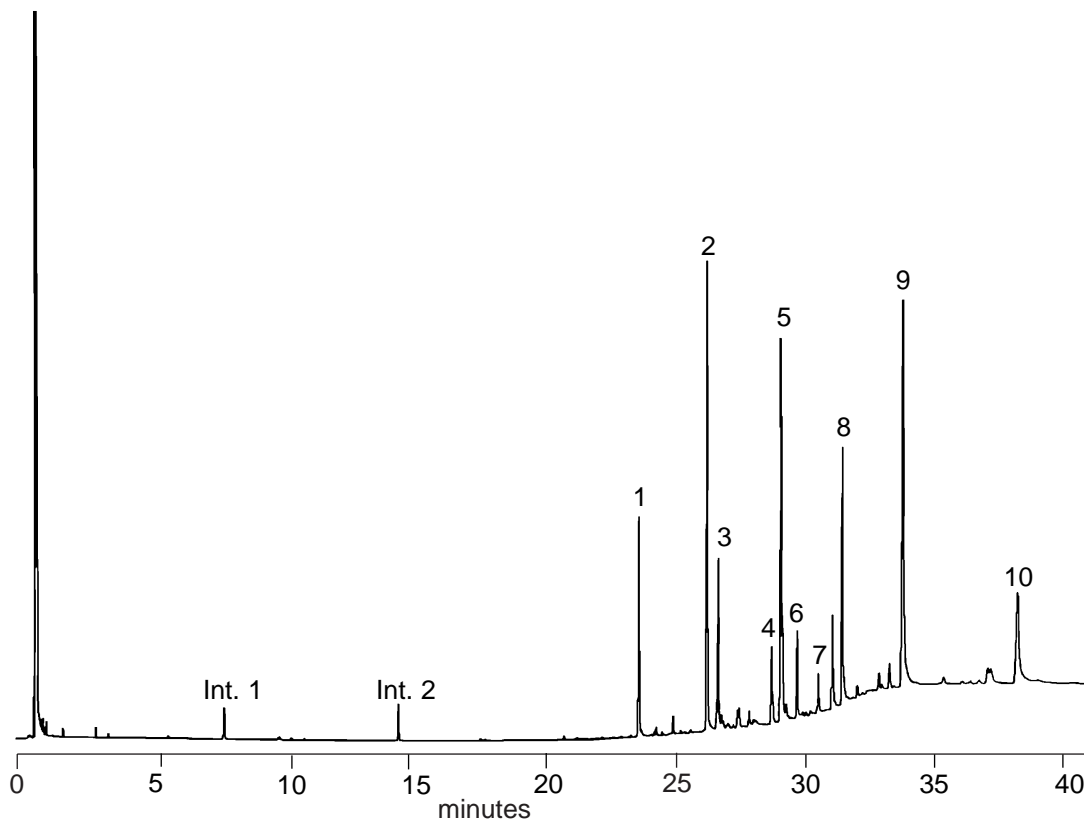

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email: support@sge.com

SGE Incorporated (USA)
email: usa@sge.com

SGE (France) Sarl
email: france@sge.com

SGE Deutschland GmbH
email: germany@sge.com

SGE Japan Inc.
email: sgejpn@itjit.ne.jp

SGE Australia Pty. Ltd.
email: sales@sge.com

SGE Europe Ltd. (UK)
email: uk@sge.com

SGE (Italia) Srl.
email: sge.italia@tin.it

SGE China Service Centre
email: sge@bj.col.com.cn

SGE India
email: sgeindia@bom4.vsnl.net.in