

PHA 08

PHARMACEUTICALS/DRUGS


Steroid Analysis (underivatised) on BPX5

STERIODS ANALYSIS (UNDERIVATISED)

Components

Part No.: 054113
Phase: BPX5, 0.25µm
Column: 25m x 0.22mm ID
Initial Temp.: 180°C
Rate: 8°C/min
Final Temp.: 350°C, 10 min
Detector: FID
Sensitivity: 32 x 10⁻¹² AFS
Injection Mode: Split
Carrier Gas: H₂, 10psi

1. Testosterone
2. Pregnenolone
3. Progesterone
4. Cholesterol

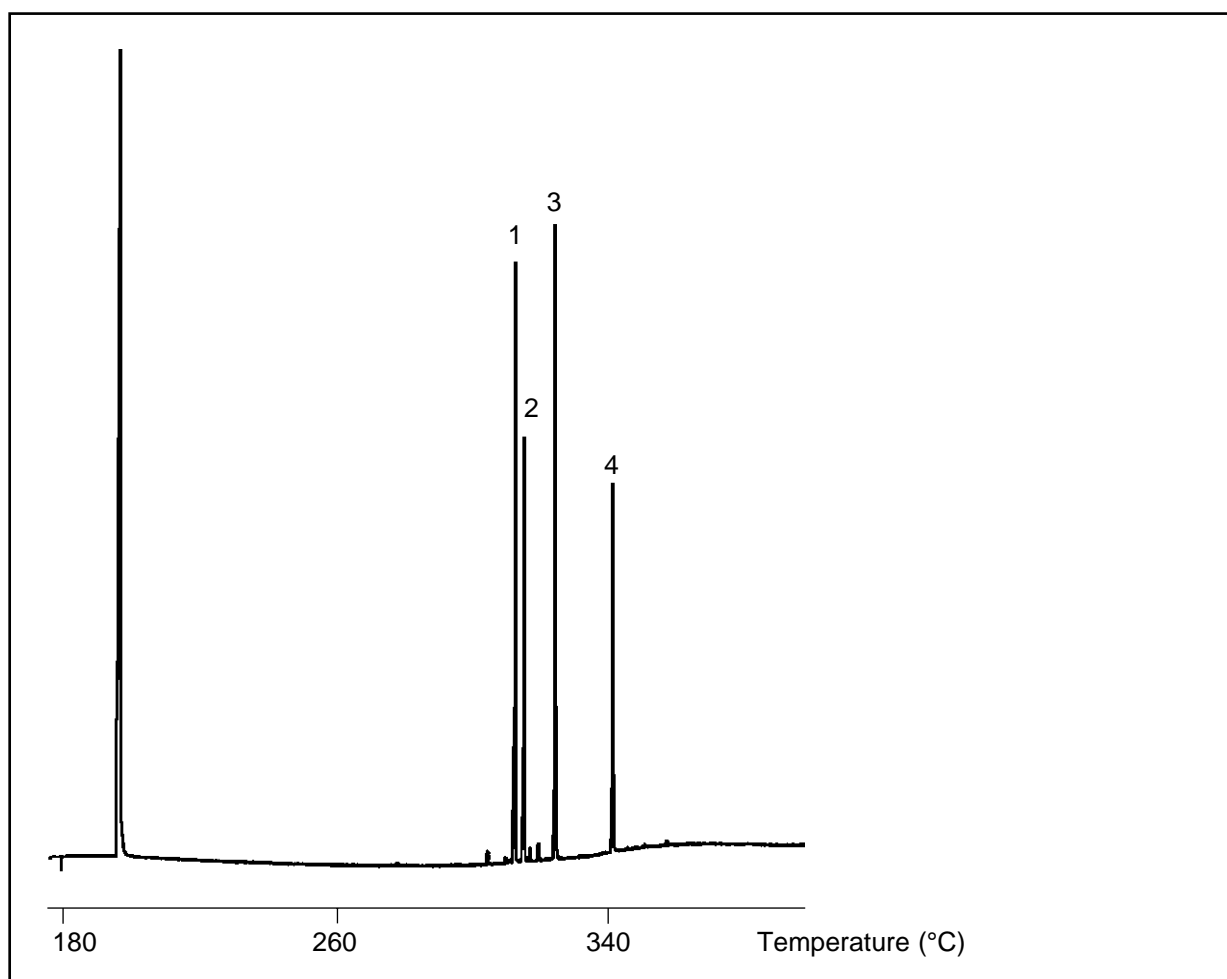

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email: support@sge.com

SGE Incorporated (USA)
email: usa@sge.com

SGE (France) Sarl
email: france@sge.com

SGE Deutschland GmbH
email: germany@sge.com

SGE Japan Inc.
email: sgejpn@itjit.ne.jp

SGE Australia Pty. Ltd.
email: sales@sge.com

SGE Europe Ltd. (UK)
email: uk@sge.com

SGE (Italia) Srl.
email: sge.italia@tin.it

SGE China Service Centre
email: sge@bj.col.com.cn

SGE India
email: sgeindia@bom4.vsnl.net.in