

PHA 10

PHARMACEUTICALS/DRUGS


## Underivatised Barbiturates on BP5

### UNDERIVATISED BARBITURATES

Part No.: 054197  
Phase: BP5, 1.0µm  
Column: 12m x 0.53mm I.D.  
Temp: 195°C  
Carrier Gas: Hydrogen  
Carrier Flow: 10mL/min  
Injection Volume: 0.1µL

### Components

1. Barbitol
2. Butabarbital
3. Amobarbital
4. Pentabarbital
5. Secobarbital
6. Hexabarbital

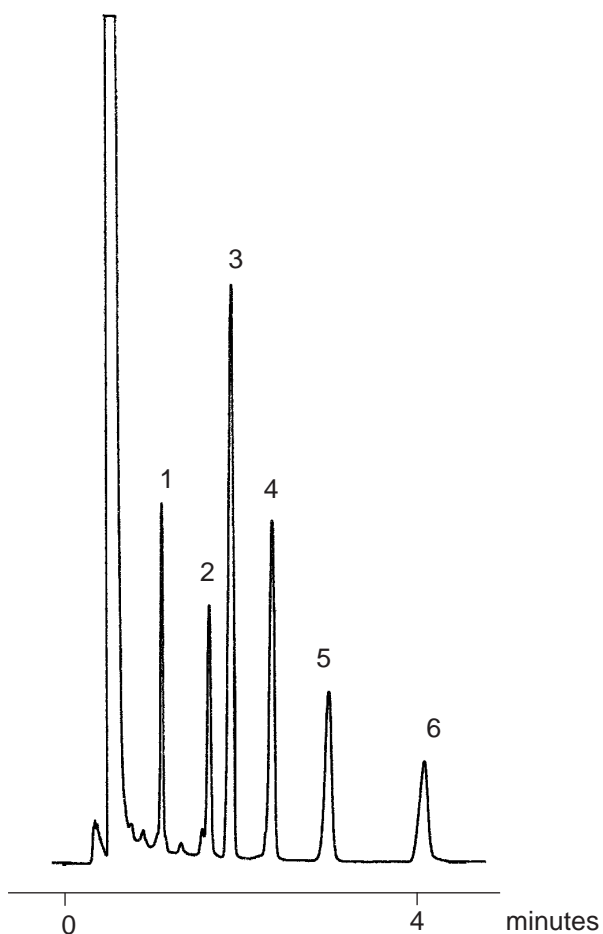

visit [www.sge.com](http://www.sge.com) or contact our technical Customer Support Team on: [techsupport@sge.com](mailto:techsupport@sge.com)

SGE International Pty. Ltd.  
email: support@sge.com

SGE Incorporated (USA)  
email: usa@sge.com

SGE (France) Sarl  
email: france@sge.com

SGE Deutschland GmbH  
email: germany@sge.com

SGE Japan Inc.  
email: sgejpn@itjit.ne.jp

SGE Australia Pty. Ltd.  
email: sales@sge.com

SGE Europe Ltd. (UK)  
email: uk@sge.com

SGE (Italia) Srl.  
email: sge.italia@tin.it

SGE China Service Centre  
email: sge@bj.col.com.cn

SGE India  
email: sgeindia@bom4.vsnl.net.in