

HPLC ANALYTICAL SOLUTIONS




LC-ACI 02
Acids

4 Carboxylic Acids on W5C18 HG

Part No.: 2066051
Packing: Wakosil-II 5C18 HG
Column: 150mm x 4.6mm
Mobile Phase: 40% Methanol, 60% 20mM H₃PO₄ adjusted to pH 3.0 with NaH₂PO₄ (v/v)
Flow: 1.0 mL/min
Temperature: Ambient
Detector: UV at 210nm

Components

1. 1,2,4,5-Benzenetetracarboxylic Acid
2. Trimesic Acid
3. Phthalic Acid
4. Benzoic Acid

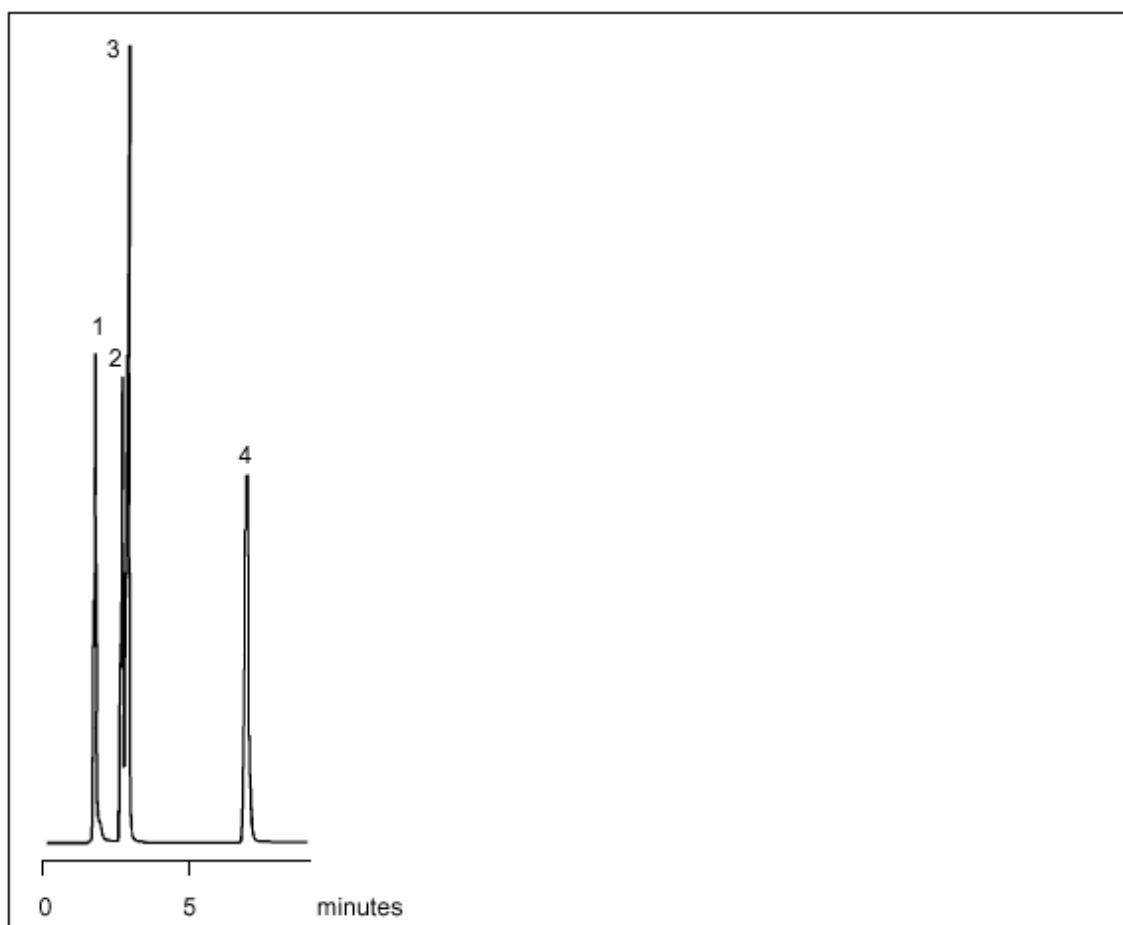

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email: support@sge.com
SGE Australia Pty. Ltd.
email: sales@sge.com

SGE Incorporated (USA)
email: usa@sge.com
SGE Europe Ltd. (UK)
email: uk@sge.com

SGE (France) Sarl
email: france@sge.com
SGE (Italia) Srl.
email: sge.italia@tin.it

SGE Deutschland GmbH
email: germany@sge.com
SGE China Service Centre
email: sge@bj.col.com.cn

SGE Japan Inc.
email: sgejpn@itjit.ne.jp
SGE India
email: sgeindia@bom4.vsnl.net.in