


**LC-ACI 05**  
**Acids**

## 10 Aromatic Organic Acids on W5C8 RS

**Part No.:** 206610  
**Packing:** Wakosil-II 5C8 RS  
**Column:** 150mm x 4.6mm  
**Mobile Phase:** 10% Acetonitrile, 90% 50mM  
KH<sub>2</sub>PO<sub>4</sub>, H<sub>3</sub>PO<sub>4</sub> (pH 4.0)  
**Flow:** 1.0 mL/min  
**Temperature:** 40°C  
**Detector:** UV at 254nm

### Components

1. 3,4-Dihydroxyphenylacetic Acid
2. Gentisic Acid
3. p-Hydroxyphenylacetic Acid
4. Salicylic Acid
5. Vanillic Acid
6. Caffeic Acid
7. p-Coumaric Acid
8. m-Coumaric Acid
9. Ferulic Acid
10. o-Coumaric Acid

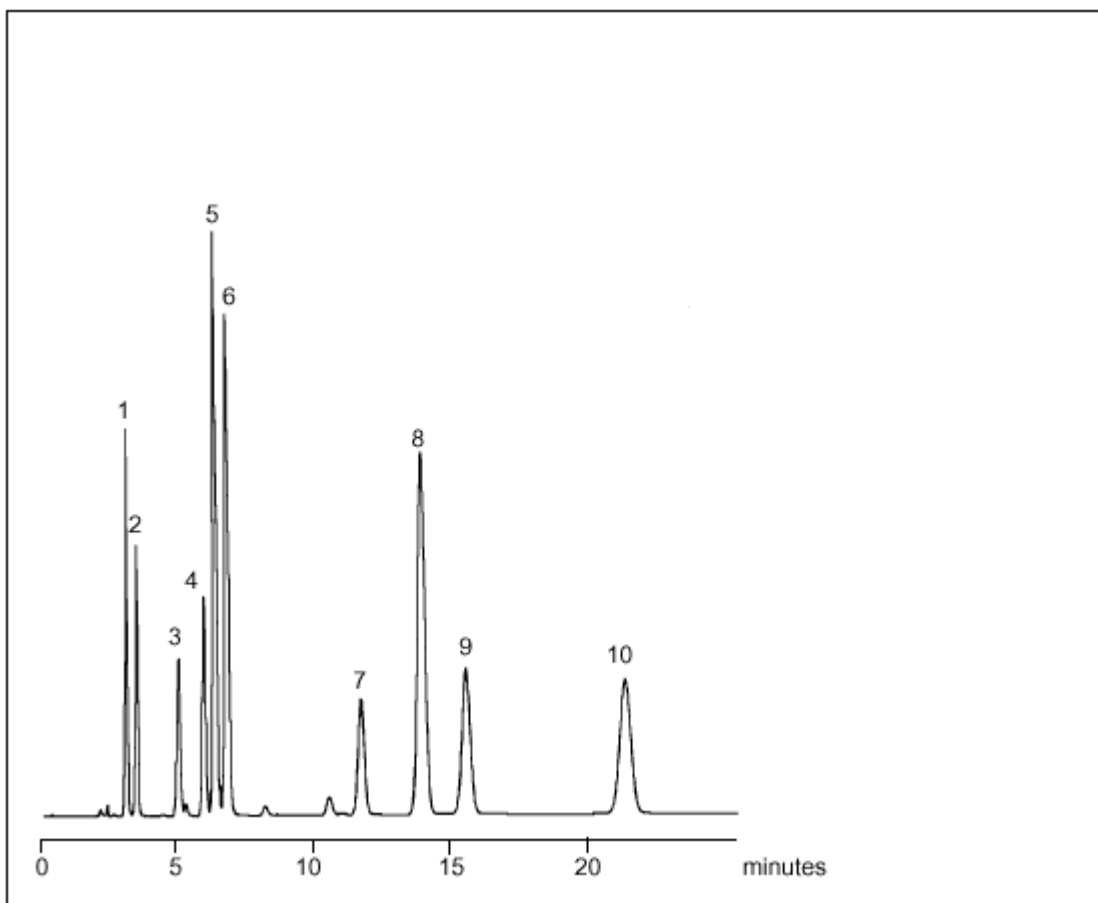

visit [www.sge.com](http://www.sge.com) or contact our technical Customer Support Team on: [techsupport@sge.com](mailto:techsupport@sge.com)

SGE International Pty. Ltd.  
email: support@sge.com

SGE Incorporated (USA)  
email: usa@sge.com

SGE (France) Sarl  
email: france@sge.com

SGE Deutschland GmbH  
email: germany@sge.com

SGE Japan Inc.  
email: sgejpn@itjit.ne.jp

SGE Australia Pty. Ltd.  
email: sales@sge.com

SGE Europe Ltd. (UK)  
email: uk@sge.com

SGE (Italia) Srl.  
email: sge.italia@tin.it

SGE China Service Centre  
email: sge@bj.col.com.cn

SGE India  
email: sgeindia@bom4.vsnl.net.in