


LC-ACI 06
Acids

10 Aromatic Organic Acids on W5C18 RS

Part No.: 206605
Packing: Wakosil-II 5C18 RS
Column: 150mm x 4.6mm
Mobile Phase: 10% Acetonitrile, 90% 50mM
KH₂PO₄, H₃PO₄ (pH 4.0)
Flow: 1.0 mL/min
Temperature: 40°C
Detector: UV at 254nm

Components

1. 3,4-Dihydroxyphenylacetic Acid
2. Gentisic Acid
3. p-Hydroxyphenylacetic Acid
4. Salicylic Acid
5. Vanillic Acid
6. Caffeic Acid
7. p-Coumaric Acid
8. m-Coumaric Acid
9. Ferulic Acid
10. o-Coumaric Acid

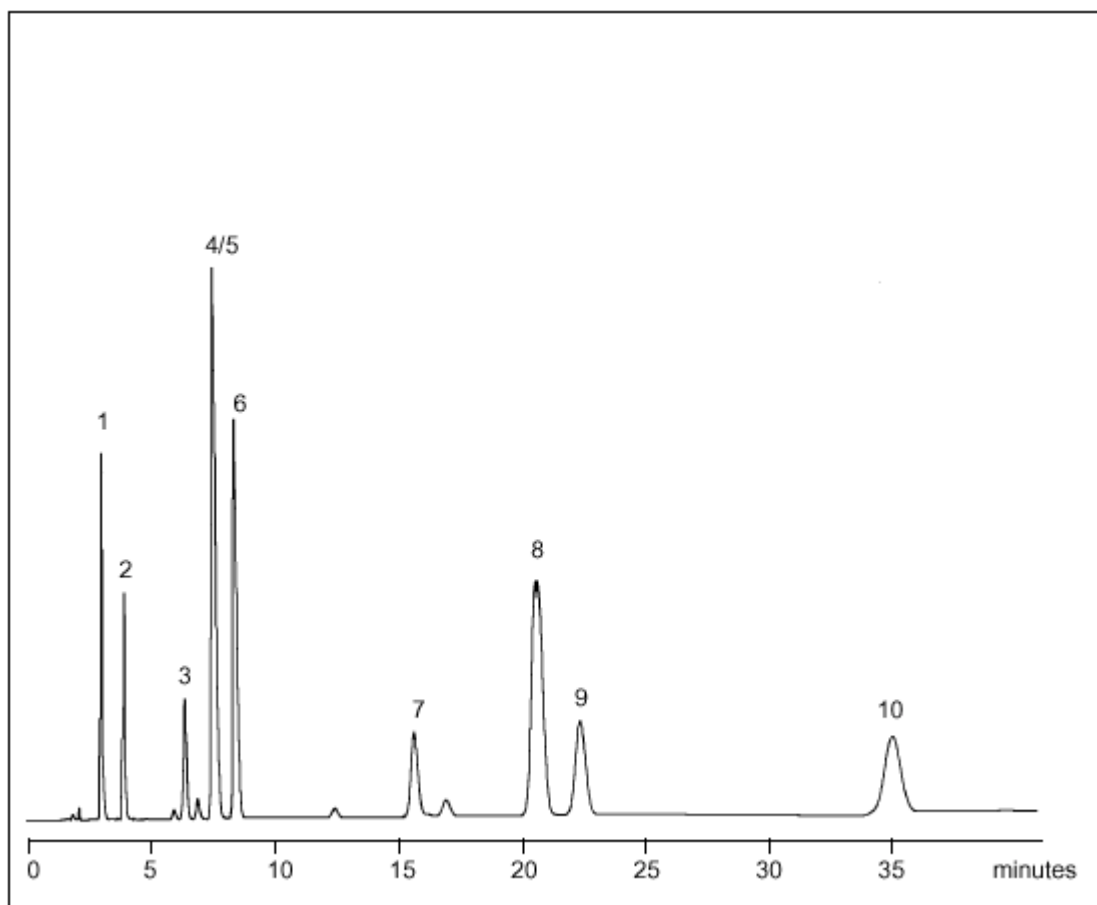

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email:support@sge.com

SGE Incorporated (USA)
email:usa@sge.com

SGE (France) Sarl
email:france@sge.com

SGE Deutschland GmbH
email:germany@sge.com

SGE Japan Inc.
email:sgejpn@itjit.ne.jp

SGE Australia Pty. Ltd.
email:sales@sge.com

SGE Europe Ltd. (UK)
email:uk@sge.com

SGE (Italia) Srl.
email:sge.italia@tin.it

SGE China Service Centre
email:sge@bj.col.com.cn

SGE India
email:sgeindia@bom4.vsnl.net.in