


LC-ARO 13
Aromatics

Pyridine on W5C18 RS

Part No.: 206605
Packing: Wakosil-II 5C18 RS
Column: 150mm x 4.6mm
Mobile Phase: 30% Acetonitrile,
70% Water (v/v)
Flow: 1.0 mL/min
Temperature: 30°C
Detector: UV at 254nm

Components
1. Pyridine
2. Phenol

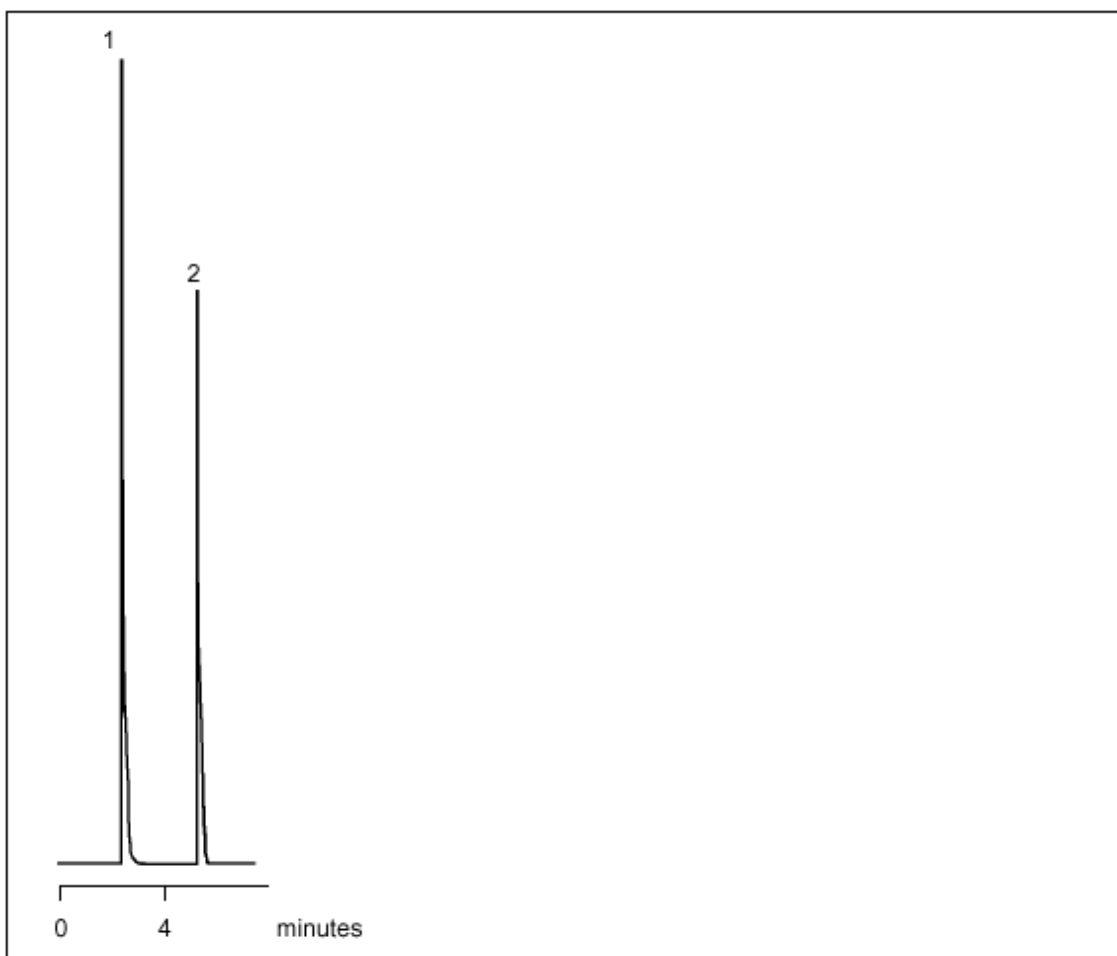

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email: support@sge.com
SGE Australia Pty. Ltd.
email: sales@sge.com

SGE Incorporated (USA)
email: usa@sge.com
SGE Europe Ltd. (UK)
email: uk@sge.com

SGE (France) Sarl
email: france@sge.com
SGE (Italia) Srl.
email: sge.italia@tin.it

SGE Deutschland GmbH
email: germany@sge.com
SGE China Service Centre
email: sge@bj.col.com.cn

SGE Japan Inc.
email: sgejpn@itjit.ne.jp
SGE India
email: sgeindia@bom4.vsnl.net.in