

LC-ARO 14 Aromatics


Aromatic Hydrocarbons on W5C18 AR

Keywords: PAH, Environmental

Part No.: 2065050
Packing: Wakosil-II 5C18 AR
Column: 250mm x 4.6mm
Mobile Phase: 85% Acetonitrile,
15% Water (v/v)
Flow: 1.0 mL/min
Temperature: Ambient
Detector: UV at 254nm, 0.16 AUFS

Components

1. Benzo(e)pyrene
2. Benzo(b)fluoranthene
3. Perylene
4. Benzo(k)fluoranthene
5. Benzo(a)pyrene

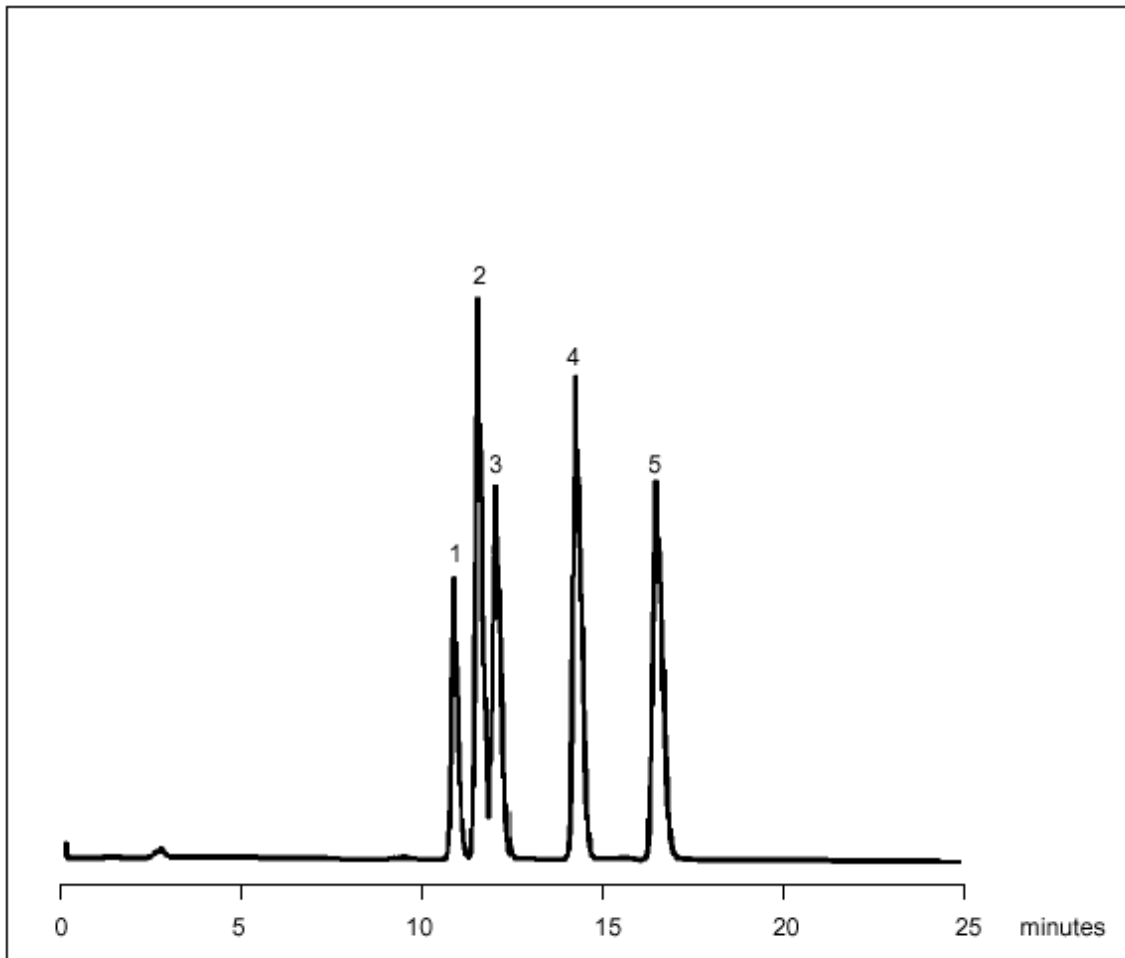

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email: support@sge.com

SGE Incorporated (USA)
email: usa@sge.com

SGE (France) Sarl
email: france@sge.com

SGE Deutschland GmbH
email: germany@sge.com

SGE Japan Inc.
email: sgejpn@itjit.ne.jp

SGE Australia Pty. Ltd.
email: sales@sge.com

SGE Europe Ltd. (UK)
email: uk@sge.com

SGE (Italia) Srl.
email: sge.italia@tin.it

SGE China Service Centre
email: sge@bj.col.com.cn

SGE India
email: sgeindia@bom4.vsnl.net.in