


LC-ARO 20
Aromatics

Pyridine Derivatives on W5C18 HG

Part No.: 2066051
Packing: Wakosil-II 5C18 HG
Column: 150mm x 4.6mm
Mobile Phase: 25% Acetonitrile, 75% 10mM K₂HPO₄ at pH 7.0 (v/v)
Flow: 1.0 mL/min
Temperature: 30°C
Detector: UV at 254nm, 0.16 AUFS

Components

1. Pyridine
2. Phenol
3. 2,4,6-Trimethylpyridine
4. 2-Propylpyridine

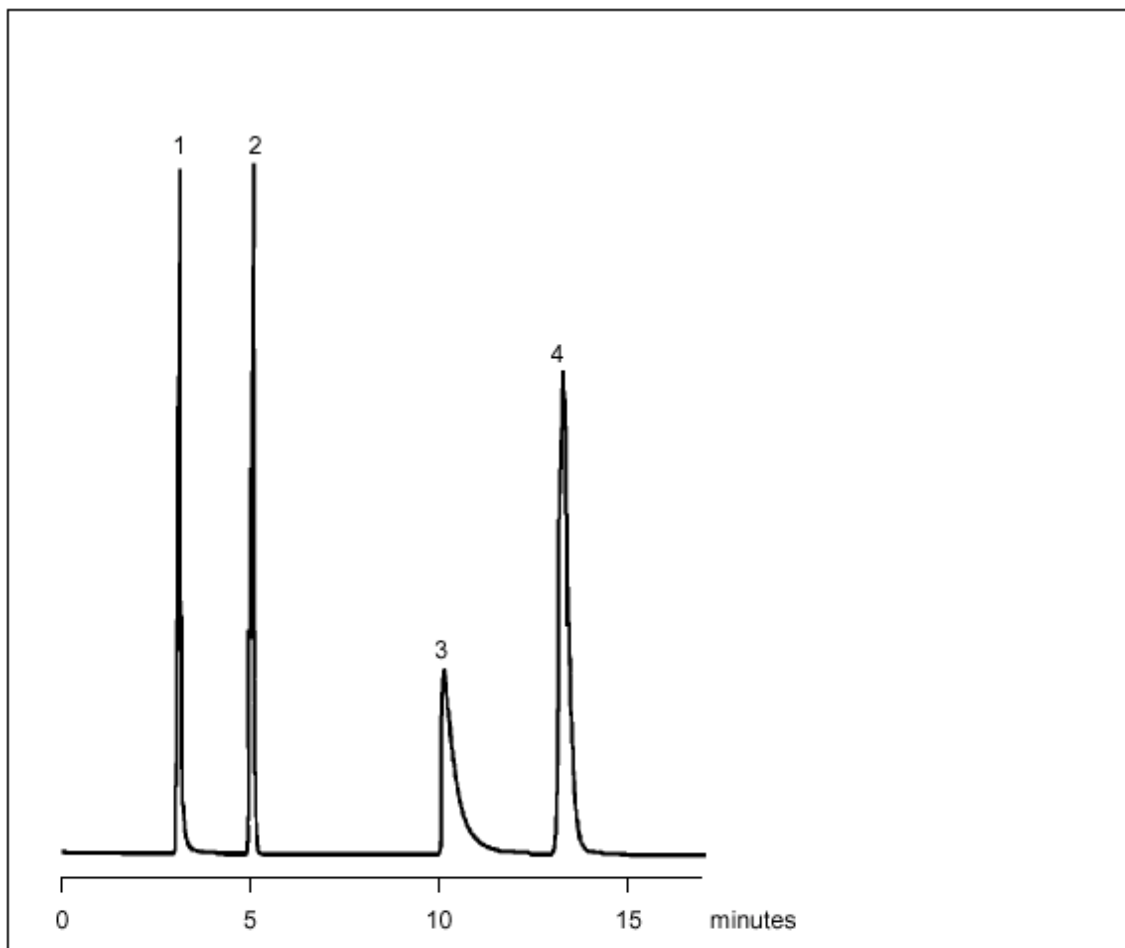

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email: support@sge.com

SGE Incorporated (USA)
email: usa@sge.com

SGE (France) Sarl
email: france@sge.com

SGE Deutschland GmbH
email: germany@sge.com

SGE Japan Inc.
email: sgejpn@itjit.ne.jp

SGE Australia Pty. Ltd.
email: sales@sge.com

SGE Europe Ltd. (UK)
email: uk@sge.com

SGE (Italia) Srl.
email: sge.italia@tin.it

SGE China Service Centre
email: sge@bj.col.com.cn

SGE India
email: sgeindia@bom4.vsnl.net.in