

LC-FOO 02


Food & Beverages

Water Soluble Vitamins on W5C18 AR

Part No.: 2066050
Packing: Wakosil-II 5C18 AR
Column: 150mm x 4.6mm
Mobile Phase: 10% Acetonitrile, 90%
(0.1% H₃PO₄, 5mM Hexane
Sulphonic Acid Na salt (v/v))
Flow: 1.0 mL/min
Temperature: 40°C
Detector: UV at 210nm

Components

1. L-(+) Ascorbic Acid
2. Nicotinic Acid
3. Nicotinamide
4. Pyridoxine Hydrochloride
5. Caffeine
6. Thiamine Hydrochloride
7. D-(+) Biotin
8. Riboflavin

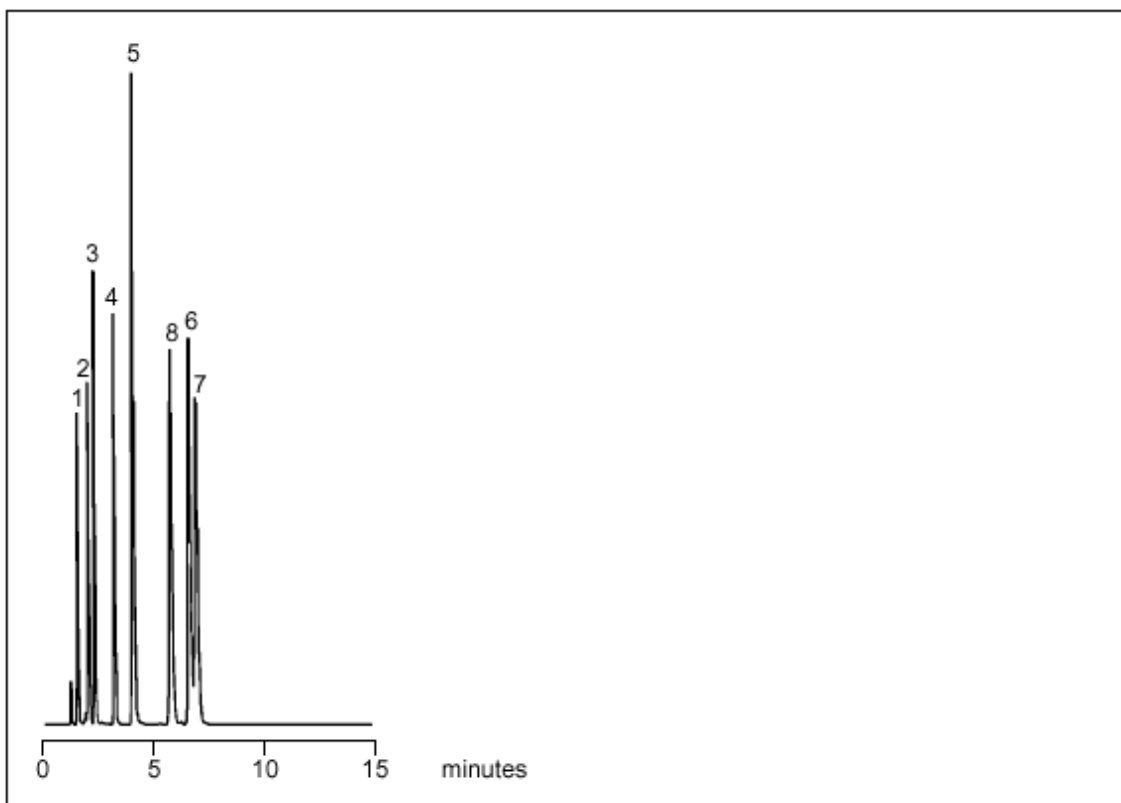

visit www.sge.com or contact our technical Customer Support Team on: techsupport@sge.com

SGE International Pty. Ltd.
email: support@sge.com

SGE Incorporated (USA)
email: usa@sge.com

SGE (France) Sarl
email: france@sge.com

SGE Deutschland GmbH
email: germany@sge.com

SGE Japan Inc.
email: sgejpn@itjit.ne.jp

SGE Australia Pty. Ltd.
email: sales@sge.com

SGE Europe Ltd. (UK)
email: uk@sge.com

SGE (Italia) Srl.
email: sge.italia@tin.it

SGE China Service Centre
email: sge@bj.col.com.cn

SGE India
email: sgeindia@bom4.vsnl.net.in