

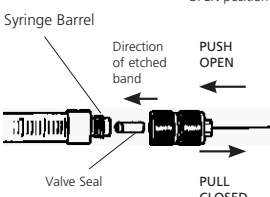
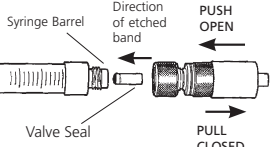
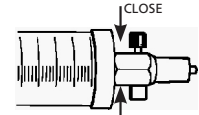
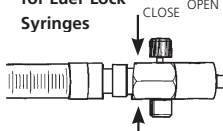
SYRINGE VALVES - GAS AND LIQUID TIGHT

Cleaning Instructions

Clean valve immediately after use by rinsing with suitable solvent (usually acetone), then dry.

Valve Information

Do not soak entire valve for long period of time.

Valve Type	To fit syringe	Valve P/N	Description	Needle Replacement P/N
<p>Push-Pull Valve with needle</p> <p>Valve shown in OPEN position</p> 	25 µL - 2.5 mL	031905 031906	V25/2.5M-0.63C V25/2.5M-0.47C	<ul style="list-style-type: none"> 038810 (0.63 cone) 038815 (0.63 sidehole dome) 038820 (0.47 cone)
<p>Push-Pull Valve with Luer Lock</p> <p>Valve shown in CLOSED position</p> 	25 µL - 2.5 mL	031907	VLL25/2.5M	Luer lock needles
<p>Push-button Valve with Luer Lock</p> <p>Valve shown in OPEN position</p> 	5 mL - 2 L	031910	VLLMA	Luer lock needles
<p>Push-button Valve for Luer Lock Syringes</p> <p>Valve shown in OPEN position</p> 	Luer tip & Luer lock 50 µL - 2 L	031915	SLLV	Luer lock needles

*Replacement Valve seal assembly (P/N 031912)

*Replacement Valve seal assembly (P/N 031912)

**CAUTION**

Particulate matter will cause the valve to become stiff or seize.

PLEASE DO NOT RETURN ANY VALVE CONTAMINATED WITH RADIOACTIVE, INFECTIOUS OR HAZARDOUS MATERIALS.

WARRANTY

SGE valves are warranted to be free of material and workmanship defects. No other warranty is expressed or implied. Any valve failing due to such defects will be replaced free of charge provided proof is given.

AUSTRALIA & PACIFIC REGION

SGE Analytical Science Pty Ltd
Toll Free: 1800 800 167
Tel: +61 (0) 3 9837 4200
Fax: +61 (0) 3 9874 5672
Email: support@sge.com

CHINA

SGE Shanghai Representative Office
Tel: +86 21 6407 9382
Fax: +86 21 6407 9386
Email: china@sge.com

MIDDLE EAST

SGE Gulf
Tel: +971 6 557 3341
Fax: +971 6 557 3541
Email: gulfsupport@sge.com

EUROPE

SGE Europe Ltd
European Head Office
Toll Free: 00800 2790 8999
Toll Free Fax: 00800 2626 2609
Tel: +44 1908 568 844
Fax: +44 1908 566 790
Tel France: +33 1 69 29 80 90
Fax France: +33 1 69 29 09 25
Tel Germany: +49 (0) 6155 / 60746 0
Fax Germany: +49 (0) 6155 / 60746 50
Email: europe@sge.com

INDIA

SGE Laboratory Accessories Pvt Ltd
Tel: +91 22 24715896
Fax: +91 22 24716592
Email: sgeindia@vsnl.com

JAPAN

SGE Japan Inc
Tel: +81 45 222 2885
Fax: +81 45 222 2887
Email: japan@sge.com

UNITED STATES OF AMERICA

SGE Incorporated
Toll Free: (800) 945 6154
Tel: +1 512 837 7190
Fax: +1 512 836 9159
Email: usa@sge.com

