


ANALYSIS OF 15 ORGANIC ACIDS ON BP20

ORGANIC ACIDS

Column Part No.: 054433

Phase: BP20, 0.25 μm
Column: 30 m x 0.32 mm ID
Initial Temp.: 70 $^{\circ}\text{C}$
Rate: 10 $^{\circ}\text{C}/\text{min}$
Final Temp.: 260 $^{\circ}\text{C}$, 5 min
Detector: FID
Injection Mode: Split
Carrier Gas: H_2 , 6 psi

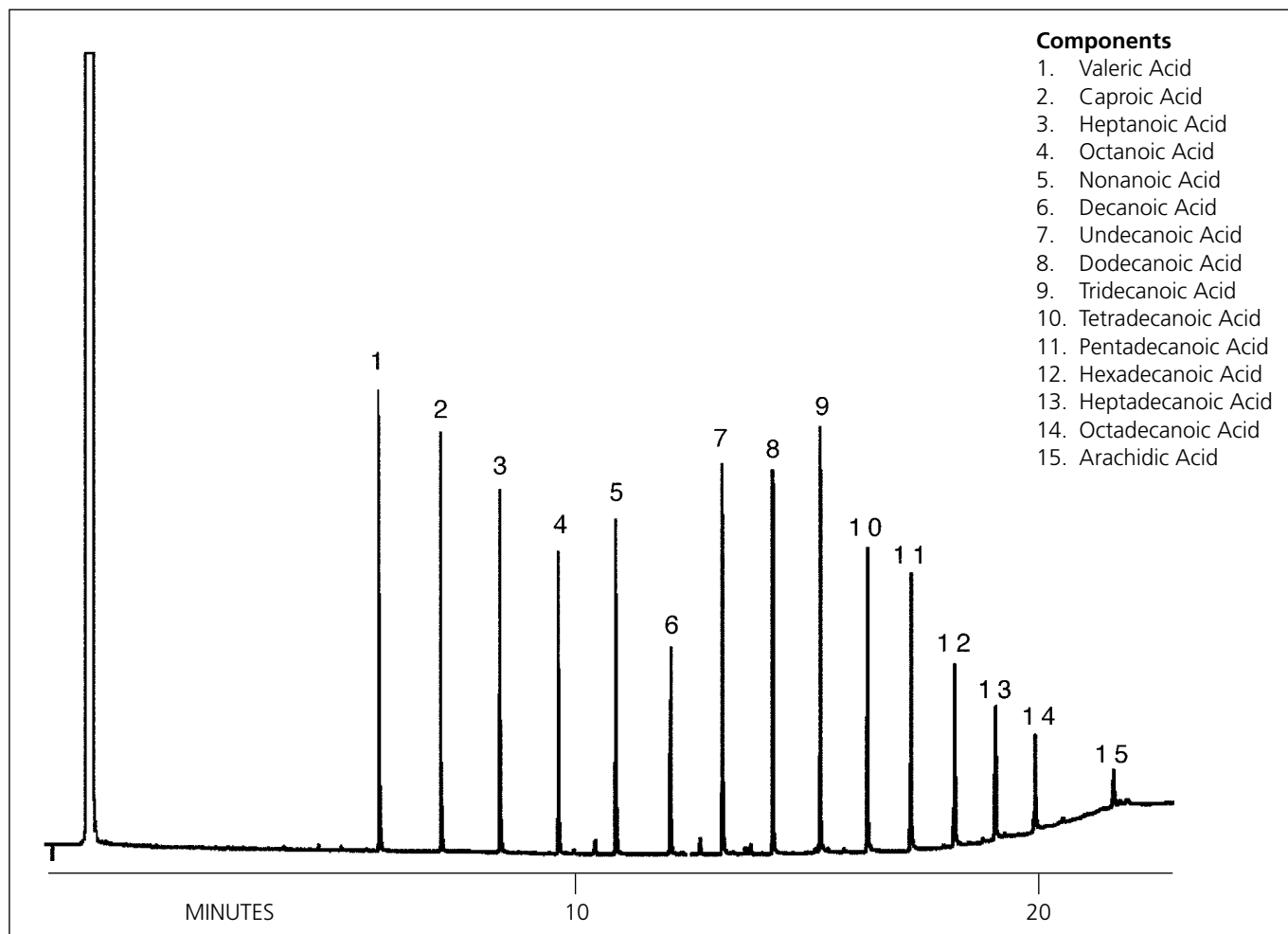

For more information contact our technical customer support team on: techsupport@sge.com