

Accessories | Gas and Vapor Sampler

The Gas and Vapor sampler (GAV-200) is ideal for remote sampling of environmental or low pressure stream gases and vapor.

Samples are collected in a 200 cm³ self-sealing, impervious bag which can be removed from the sampling canister for further sampling. The filled bags can be stored until discharge of the contents is required.

Features and Benefits

- Easy to use.
- Ideal for field environmental analysis.
- Resealable/removable sample bags.
- 200 cm³ sample bag capacity.

Recommended Applications


Spot collection of environmental gas samples or low pressure gas streams.



Product Specifications

Kit Contains:

- Canister,
- Bulb,
- 10 bags,
- Marker Pen, and
- Instructions.

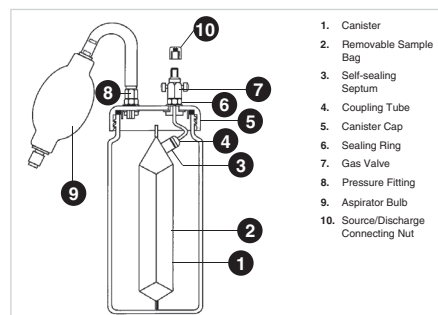

The GAV-200 gas sampling bags consist of a polyethylene/aluminum foil/polyethylene laminate. The sample bags are suitable for gas pressures differentials up to 5psi. The canister is polycarbonate with a polypropylene cap. The maximum canister pressure is 50 psi (350 kPa).

For more information about this product visit www.sge.com or contact techsupport@sge.com.